## Stage 1


- 1. Weigh out 0.500 g aminothiazole
- 2. Dissolve with 10 cm<sup>3</sup> THF in a vessel
- 3. Cool it using an ice bath
- 4. While there, add 0.15 cm<sup>3</sup> lodomethene (micro)dropwise; very slowly
- 5. Put in a stopper (prevent evaporation) and place in fridge.

AMINOTHIAZOLE

0

## @

①

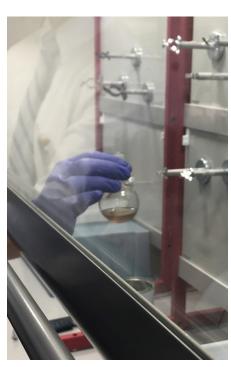





- 1. Cross fingers
- 2. Mix with a -pyrodine

## Stage 3

- 1. Move to separating funnel
- 2. Add salt water and mix well; allow to separate
- 3. Tap off the aqueous solution and the organic solution into separate vessels

